



You made my night.



The ResMed AirFit F40

A minimalist mask.
Maximized for everyone.

AirFit™
F40

Introducing the ResMed AirFit™ F40

This mask builds on what you—and your patients—love about ResMed masks, adding a fully flexible frame with new soft silicone cushions and fewer adjustment points for easy fitting and comfort. An UltraCompact under-the-nose design allows flexible bedtime routines and sleeping positions. And, it's easy to snap on and snap off, even in the middle of the night.

98%
Fit range¹

New AdaptiSeal™ full face cushion

Soft silicone adapts to the face, designed to maintain seal.

Quick-release elbow

Makes it easy to connect and disconnect to therapy.

Short tube

Reduces drag and moves with the patient.

QuietAir™ diffused venting

Integrated into the cushion to provide a more peaceful bedroom environment.

Headgear adjustment points 1 2 3 4

Two upper and two lower adjustment points simplifies fitting and use.

Flexible mask frame

Plastic frame covered by a soft nylon fabric padding. The frame sits below the eyes, designed for comfort and to reduce facial marks.

Magnetic clips*

Easily snaps on and off, even in the middle of the night.



*The mask contains magnets that may interfere with certain implants or medical devices. Please refer to the User Guide for complete labeling information, including magnet contraindications and warnings.



Complete mask system

- 64600
- 64601
- 64602
- 64603
- 64604



AirFit™ F40

Frame system 64607, 64608, 64609

AirFit F40

Complete mask system

(includes mask frame, full face cushion and headgear)

Product name	Model number
AirFit F40 System – SW/STD AMER	64600
AirFit F40 System – M/STD AMER	64601
AirFit F40 System – L/STD AMER	64602
AirFit F40 System – M/SML AMER	64603
AirFit F40 System – M/LGE AMER	64604

Frame system

(includes mask frame and full face cushion)


Product name	Model number
AirFit F40 Frame System – SW	64607
AirFit F40 Frame System – M	64608
AirFit F40 Frame System – L	64609

Full face cushion

Product name	Model number
AirFit F40 Cushion Small Wide	64627
AirFit F40 Cushion Medium	64628
AirFit F40 Cushion Large	64629

Headgear

Product name	Model number
AirFit F40 Headgear SML	64633
AirFit F40 Headgear STD	64634
AirFit F40 Headgear LGE	64635

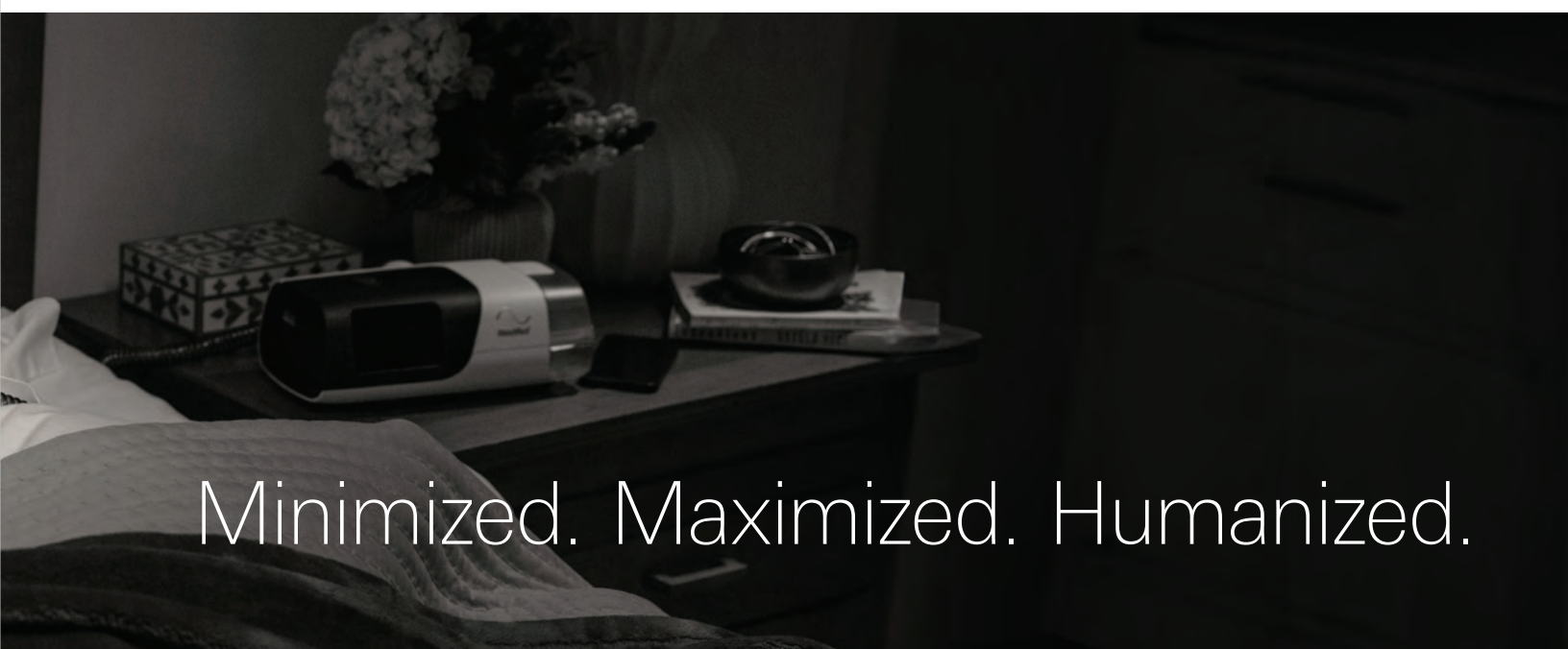

Minimized. Maximized. Humanized.

How to fit the AirFit F40



1
Unclip the lower straps away from the frame and place cushion under the nose to fit against the face.
With the ResMed logo facing up, pull headgear over head.



2
Bring the lower straps behind the ears and reattach to the frame with the magnetic clips*
Feed longer hair through the back of the frame, if needed.



3
Adjust the upper straps and then the lower straps for a snug fit.



Scan the QR code for full video instructions.

Experience the ResMed AirFit F40 today.

The ultra-compact, full face mask that combines high seal performance with comfort in a flexible, easy-to-fit design. And, it's brought to you by ResMed, the sleep apnea brand most trusted by sleep physicians.²

The ResMed AirFit F40. A minimalist mask, optimized for everyone.



Scan the QR code or contact your ResMed sales representative.



*This mask contains magnets that may interfere with certain implants or medical devices. Please refer to the User Guide for complete labeling information, including magnet contraindications and warnings.

¹ Two ResMed external clinical studies of experienced PAP patients. The first study included users with ≥ 6 months' therapy who used ResMed's AirFit F40 and AirFit F20 full face masks in randomized order for up to seven consecutive nights each at home. Study conducted in Australia, June 6 – 27, 2023, n=28 (AirFit F40) and 29 (AirFit F20). Facial measurements of 29 participants were taken. The second study included current users with ≥ 12 months' PAP therapy who met 90-day CMS compliance for PAP therapy in the three months prior to the study. Participants used ResMed's AirFit F40 and Fisher & Paykel's Evora Full in a randomized order with their own for up to seven consecutive nights each at home. Study conducted in the U.S., October 4 – November 20, 2023, n=57. Facial measurements of 57 participants were taken.

² ResMed study of board certified sleep physicians in the U.S. Study conducted August 12 – September 6, 2022. n=200